

Capp *Rondo*

4-PLACE PLATE SHAKER

Capp *Rondo* 4-Place Plate Shaker is the ideal choice for thorough, simultaneous mixing in up to 4 individual microplates of different types. The 3mm orbital motion movement guarantees the ultimate mixing performance and the brushless DC motor ensures long and maintenance free operation, making **Capp *Rondo* 4-Place Plate Shaker** the ideal device for extended runs.

Capp *Rondo* 4-Place Plate Shaker operates with the speed range from 200 to 1200 RPM. Its small and compact design requires minimum benchtop space, and the digital display gives the user an overview of the speed and time settings.


Capp *Rondo* 4-Place Plate Shaker can operate from 1 to 999 minutes, or without any time limits - in an infinite mode.

The rubber o-rings ensure an ultimate grip for all 4 plates. At the same time, loading and unloading is done effortlessly. Once loaded, all 4 plates can be easily accessed.

The pulse mode of **Capp *Rondo* 4-Place Plate Shaker** enables continuous mixing in clockwise and counterclockwise directions, changing in a time interval from 30 to 90 seconds. Thanks to the last run memory feature, the user can comfortably continue their work in the same mode as previously.

Capp *Rondo* 4-Place Plate Shaker is compatible with all commonly used skirted microplates and extremely easy to clean and maintain.



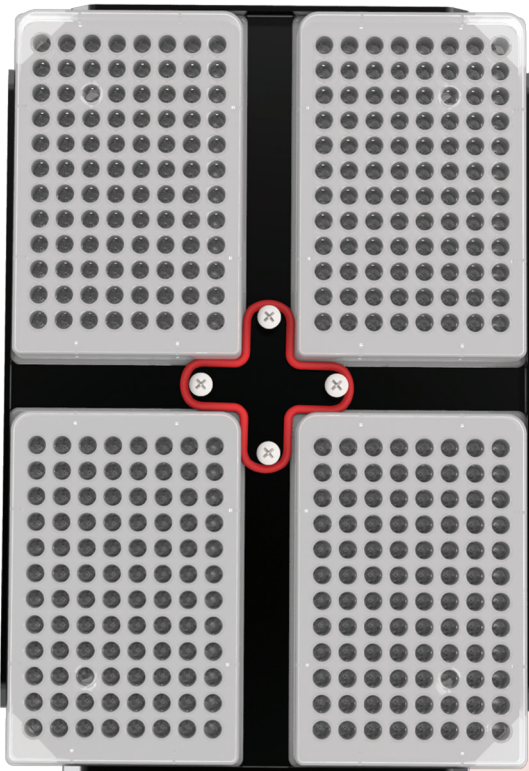

Specifications:	Capp <i>Rondo</i> 4-Place Plate Shaker
Motor type	BLDC motor
Orbital Motion	3mm
Number of plates	Up to 4
Variable Speed	200 rpm to 1200 rpm
Run Time	1 to 999 mins and infinity mode
Pulse Mode	Yes, Programmable (30 to 90 seconds)
Weight	5.4kg
Maximum Volume	4 x 384 x 60µL
Inbuilt Counterbalance	Yes
Dimensions (L x W x H)	343 x 88 x 195mm
Digital Display	Yes
Input Voltage	24V
Attachments	Acc. All standard Skirted Microplates

Ordering Information

Cat No.	Description
CRP-412X	CappRondo 4-Place Plate Shaker 1200 rpm

